



These are the standard textiles and colours. See more information on our website: [www.vela-medical.com](http://www.vela-medical.com)

## Compound



BLACK - 825

### Compound

**Cleaning:** Wipe regularly with warm water and neutral detergent. Use a clean white cloth or a sponge.

**Disinfection:** Use detergents containing up to 70-85% ethanol or chlorine disinfection containing active chlorines at 1000-1200 ppm. After cleaning, rinse repeatedly with cold, clean water.



ANTHRACITE GREY - 828



BLUE - 823

## Illusion



BLACK - 58

### Illusion

**Cleaning:** Wipe regularly with warm water and neutral detergent. Use a clean white cloth or a sponge.

**Disinfection:** Use detergents containing up to 75% ethanol or chlorine disinfection containing active chlorines up to 5000 ppm. After cleaning, rinse with cold, clean water for longer durability.




GREY - 59



TILE - 60

## Dolce

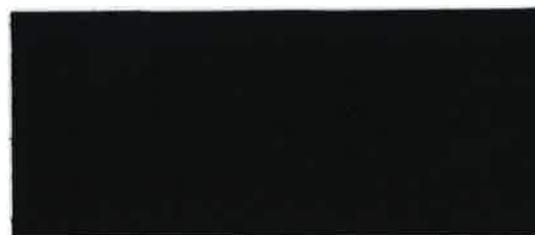

BLACK - 771

### Dolce

**Cleaning:** Wipe with a solution of 10% mild soap and warm water using a soft cloth. Rinse with a clean, damp cloth.

**Disinfection:** Use any of the following solutions: ethanol, water and bleach (5:5:2) or bleach and water (1:5) or ethanol and water (1:1).



METEOR - 777



PEACOCK - 770